

Analysis of Perchlorate in Food Samples of Plant Origin Applying the QuPPE-Method and LC-MS/MS

Julia Hepperle, Anne Wolheim, Diana Kolberg, Cristin Wildgrube, Ingrid Kaufmann-Horlacher, Michelangelo Anastassiades and Ellen Scherbaum

Chemisches und Veterinäruntersuchungsamt Stuttgart

This study reports about how to integrate the environmental contaminant perchlorate into the QuPPE multiresidue method for the analysis of polar pesticides. Perchlorate, along with a number of polar pesticides, was extracted with acidified methanol. The extracts were filtered and subjected to liquid chromatography-triple quadrupole mass spectrometry (LC-MS/MS). ¹⁸O-labeled perchlorate was used as internal standard. The assay was validated for apple and barley matrix at 0.01 and 0.1 mg/kg, showing average recoveries between 93 % and 108 % and relative standard deviations (RSDs) between 2.2 % and 5.0 %. Up to now, we have analyzed 350 conventionally and 78 organically grown food samples (vegetables, fruits, cereals, legumes, mushrooms and processed foods) from the local market for perchlorate. Perchlorate could be detected in 66 (19 %) of the conventional samples and in two (2.6 %) of the organic samples; in six samples the provisional maximum tolerable daily intake (PMTDI) of 0.01 mg/kg bodyweight was exceeded.

Introduction

In 2012, the environmental contaminant perchlorate was frequently detected in conventionally and organically grown fruit and vegetable samples from the local market. Perchlorates are salts of perchloric acid (HClO₄) arising from natural as well as from anthropogenic sources. The perchlorate anion is a stable but strong oxidant as the chlorine atom possesses an oxidation state of +7. Most perchlorates are soluble in water and circulate within the hydrologic cycle [1]. They naturally occur ubiquitous in atmospheric water and soil compartments in low concentrations [2]. However, perchlorate contamination in food cannot result from natural sources as the appearing environmental concentrations are too low. Aside from natural occurrence, perchlorates are anthropogenically spread, since they are widely used in the metalworking industry, in

paper upgrading and, as oxidants, for propellants in pyrotechnic materials like rockets, fireworks and explosives [1, 2]. In diagnostic radiology, perchlorates are medically applied in iodine contrast agents keeping the exposure to high iodine concentrations at a low level by blocking the thyroid [2, 3]. Another source of perchlorate contamination are fertilizers or sewage sludge which are both widely applied in agriculture [2]. Perchlorates also come along with chlorate, a substance which has formerly been used as a herbicide in plant protection products. However, pesticides containing chlorate as an active substance have been banned within the EU since May 2010 [4]. Like chlorate, perchlorate can act as a herbicide [1], but is not approved as a pesticide within the European Union [5] and there is no maximum residue level listed in Regulation (EC) 396/2005 [6].

Perchlorate inhibits the production of thyroxine (T₄) and triiodothyronine (T₃) in the thyroid by blocking the uptake of iodide from the blood into the thyroid [7, 8]. This effect is reversible as the perchlorate acts as a competitive inhibitor. The Joint FAO/WHO Expert Committee on Food Additives proposed a provisional maximum tolerable daily intake (PMTDI) of 0.01 mg/kg bodyweight for perchlorate [9], whereas the National Academy of Sciences (NAS) determined a Reference Dose (RfD) of 0.007 mg/kg bodyweight [10]. Both values are based on the same human toxicity study on perchlorate, in which adult human beings were exposed to perchlorate from drinking water for 14 days [8]. The study resulted in a NOEL (No Observed Effect Level) of 0.007 mg/kg bodyweight per day for perchlorate, which means that from this concentration or above, perchlorate inhibits the uptake of iodide into the thyroid.

Ion chromatography, mostly coupled to conductivity detection has been the method of choice to determine perchlorate [3, 11–13].

Other possibilities to analyze the target anion are photometry [14], potentiometry [15, 16], electrophoresis [17, 18], atomic absorption spectrometry [19, 20], and Raman spectrometry [21]. In the last few years, scientist tended to develop methods including ion chromatography or liquid chromatography coupled to mass spectrometry detection [22–24]. Up to now, the predominant subject of perchlorate analysis has been drinking water and soil [17, 18, 24–28]. To our knowledge, fewer investigations have been made in determining perchlorate in food samples [12, 29–31].

In our study, we analyzed samples of plant origin applying the QuPPE method following liquid chromatography and triple quadrupole mass spectrometry. The QuPPE-method (short for “Quick Method for the Analysis of Residues of numerous Highly Polar Pesticides in Foods of Plant Origin involving Simultaneous Extraction with Methanol and LC-MS/MS Determination”) was developed by the European Reference Laboratory for pesticide residues requiring Single Residue Methods (EURL-SRM). It includes a simple extraction with methanol followed by a filtration step. The extracts are directly analyzed via LC-MS/MS without any further cleanup. The method was validated for apple (high in water) and barley (dry) matrix and was applied to monitor conventionally and organically grown real samples.

Experimental

Chemicals and standards

Acetonitrile and methanol of gradient grade were purchased from Merck (Darmstadt, Germany). Water for liquid chromatography was deionized in the laboratory using a Direct-Q 3 UV Ultrapure Water Purification System (Billerica, MA, USA). Dry ice for sample comminution was delivered in blocks by a local provider, stored at -80 °C until use and was tested not to contain perchlorate at relevant levels. Ammonium formate was purchased from Sigma-Aldrich (Steinheim,

Germany). The analytical pesticide standard sodium perchlorate (purity > 98 %) was obtained from Sigma-Aldrich (Steinheim, Germany), the internal standard 18O-labeled perchlorate (heavy oxygen > 90 %) was purchased from Cambridge Isotope Laboratories Inc. (Andover, USA). For perchlorate, a stock solution of 1 mg/mL, considering standard purity and the sodium part, was prepared in methanol. The stock solution of the internal standard perchlorate 18O (100 µg/mL in water) was diluted to 10 µg/mL and 1 µg/mL (working solution) with methanol. All solutions were stored in the dark at 4 °C.

Samples and commodities

Different crops of plant origin were sampled from the local market. Representative parts of the fresh samples were pre-chopped and deep-frozen at -18 °C overnight after arriving in the laboratory. The pre-chopped samples were subsequently milled with dry ice and stored at -18 °C until extraction.

Dry samples (e.g. barley) were ground thoroughly without adding water and were stored at room temperature until extraction.

Apparatus

Frozen samples were chopped with dry ice using a Prime Cut UM5 universal machine from Stephan Machinery GmbH (Hameln, Germany). Dry samples were ground at room temperature using a Grindomix GM 200 knife mill by Retsch (Haan, Germany). The automatic shaking machine Geno Grinder 2010 (SPEX Sample Prep, Metuchen, USA) was used for automated extraction. The centrifuge Rotanta 460 by Hettich (Tuttlingen, Germany), appropriate to the centrifuge tubes employed in the procedure and capable of achieving 4000×g, was used. Electronic pipettes applicable for volumes of 10–100 µL and 200–1000 µL, respectively, and manual pipettes applicable for volumes of 1–10 mL were from Eppendorf (Hamburg, Germany). Analytical balances capable of weighing to 0.1 mg or to 0.01 g were from Mettler-Toledo (Greifensee, Switzerland). The Fortuna Optifix 10 mL Universal Dispenser was used to add solvent to the samples. The automatic dispenser system Opus (20–50 mL; Ex 20 °C; Hirschmann Laborgeräte, Eberstadt, Germany) was used for dilution of the working solutions.

50-mL PP (114×28 mm) single-use tubes with screw caps for the sample extraction were from Sarstedt (Nümbrecht, Germany). 1.5-mL LC autosampler vials and 20-mL screw-cap vials were from Ziemer GmbH (Mannheim, Germany). 6 mL single-use syringes from Henke Sass Wolf (Tuttlingen,

Germany) and disposable polyester syringe filters (0.45 µm pore size, 15 mm diameter) from Machery-Nagel (Düren, Germany) were used to filter the QuPPe extracts.

An Agilent 1200 HPLC system (Agilent, Waldbronn, Germany), consisting of a binary pump, a standard autosampler and a column oven, combined with the mass spectrometer API 5500 Q-Trap (ABSCIEX, Darmstadt, Germany), run in ESI negative mode, was used for the analysis of the QuPPe extracts. The samples were injected on a Hypercarb 2.1 x 100 mm column with 5 µm particle size (Thermo Fisher Scientific, Waltham, USA), equipped with a Hypercarb Guard (2.1 x 10 mm, 5 µm particle size) and a Supelco column saver 2 µm-precursor filter (Sigma-Aldrich, Steinheim, Germany).

QuPPe extraction

For the QuPPe method 10 0.1 g (in the case of cereals 5 0.1 g plus 10 mL water) of the prepared frozen samples were weighed into 50-mL PP tubes. 10 mL of acidified methanol followed by 100 µL of the internal standard working solution (10 µg/mL) were added to the samples. The tube was closed and shaken by a mechanical shaker for one minute. The tube was centrifuged for 5 min at 4000 rpm and the upper methanol phase was filtered through a syringe filter (0.45 µm) into sealable storage vessels. An aliquot of the extract was transferred into autosampler vials and was subsequently diluted with acidified methanol by a factor of 1:10 (depending on the matrix, see "Results and Discussion" below). The extracts were applied to LC-MS/MS measurement.

LC-MS/MS analysis

Mobile phase A consisted of 1 % acetic acid in purified water containing 5 % of methanol, mobile phase B was 1 % acetic acid in methanol. A gradient program was applied starting at a flow rate of 0.2 mL/min and 100 % of mobile phase A (0 % B) at injection time, gradually changing to 70 % A (30 % B) over 10 minutes. Then, the flow was set to 0.4 mL/min within 1 minute and the conditions were kept for another 7 minutes. Mobile Phase A was then shifted to 10 % (90 % B) within 1 minute and this composition was kept for 3 minutes. Then, mobile phase composition and the flow were shifted back to starting conditions within 0.1 minutes and kept for approximately 8 minutes in order to reequilibrate the column. Total run time was 30 minutes; the injection volume was 5 µL. The column temperature was set to 40 °C. ESI interface in negative mode was applied, the MS/MS detection was performed in the multiple reaction monitoring (MRM) mode.

The ion spray voltage was -4500 V, the curtain gas was nitrogen at 30 psi, the nebulizer gas and the turbo gas were synthetic air, both at 60 psi, the source temperature was 420 °C and the collision gas was set to "medium".

Two precursor-product ion transitions were monitored for perchlorate in the multiple reaction monitoring mode (MRM). The first transition (m/z 99/83) was used for quantification, while the second transition (m/z 101/85) served for confirmation. For the internal standard perchlorate 18O m/z 107/89 was monitored. Insight to additional method details is given at the web page of the European Reference Laboratory for pesticide residues requiring Single Residue Methods (EURL-SRM) [32].

Method validation

The method was validated at levels of 0.01 and 0.1 mg/kg for apple matrix (commodity with high water content) and barley matrix (dry commodity). Blank samples were spiked with perchlorate at 0.01 and 0.1 mg/kg (n=5) following QuPPe extraction and LC-MS/MS analysis. All blank samples were proved not to contain perchlorate in relevant amounts right before validation. We applied matrix matched calibration standards at concentrations representing 60 and 120 % of the validation level for each matrix. Quantification was carried out via the ratio between analyte peak area and ISTD peak area for both, calibration standards and the extracts of the recovery samples. Results of ongoing validations can be obtained via the Datapool of the EU Reference Laboratories for Residues of Pesticides [33].

Results and Discussion

In LC-MS/MS analysis we applied a porous graphitic carbon column for chromatographic separation. By using this column it was possible to avoid laborious washing procedures coming along with ion chromatography. The Thermo Hypercarb column consists of 100 % porous graphitic carbon with no micropores or chemically bonded phases. The carbon atoms are hexagonally arranged and form flat sheets. Highly polar analytes or charged compounds induce dipoles at the graphic surface and interact with the induced partial charges. Perchlorate showed a narrow and well-shaped peak at a retention time of 12.20 minutes. Figure 1 and 2 show typical chromatograms of perchlorate (0.01 mg/kg) and the 18O-labeled internal standard (0.05 µg/mL) in a QuPPe extract from apples.

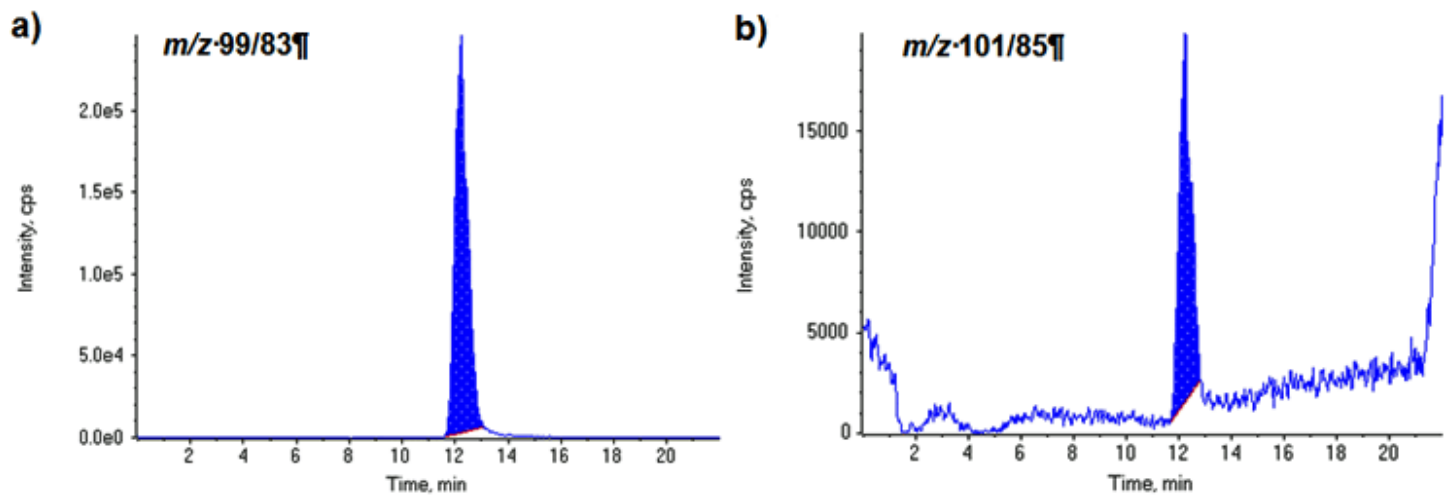


Fig. 1: LC-MS/MS chromatogram of perchlorate in apple matrix (0.01 mg/kg) a) m/z 99/83 b) m/z 101/85

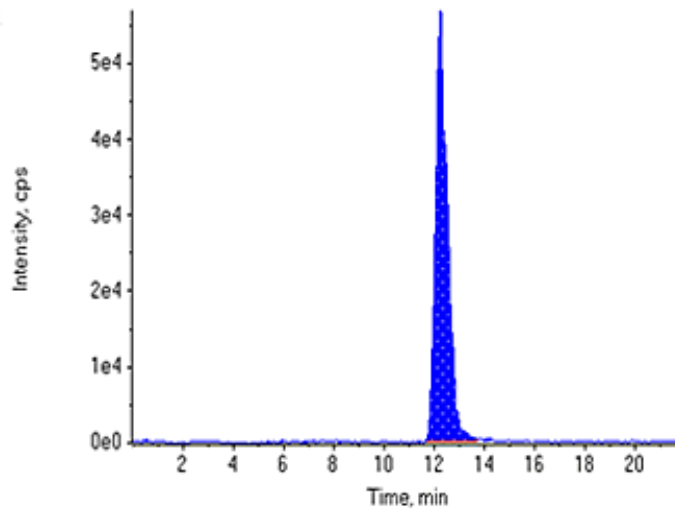


Fig. 2: LC-MS/MS chromatogram of ^{18}O -labeled perchlorate m/z 107/89 (internal standard) in apple matrix (0.05 $\mu\text{g/mL}$)

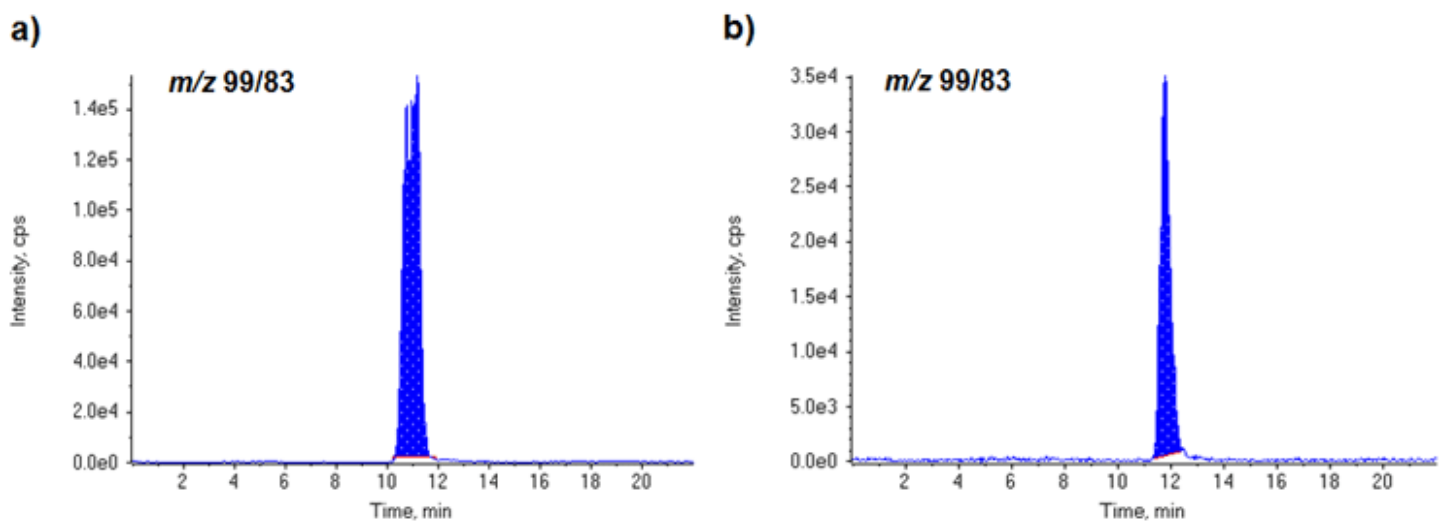


Fig. 3: LC-MS/MS chromatogram of perchlorate m/z 99/83 in kale a) without dilution b) 1:10-diluted extract

Being a strong solvent, methanol tends to co-extract matrix ingredients that may cause matrix effects in instrumental analysis. Diluting the QuPPE extracts results in less co-extracted matrix ingredients per injection and matrix effects in LC-MS/MS analysis can be avoided or reduced thereby [34]. In our case, matrix effects should be eliminated by the use of the ^{18}O -labeled internal standard. However, we observed unshaped and uneven peaks for perchlorate in certain matrices while method development. We noticed improved peak shapes in those matrices when QuPPE extracts were diluted with acidified methanol. Figure 3 shows chromatograms of a diluted and an undiluted QuPPE extract from kale containing perchlorate.

The limit of quantification (LOQ) was estimated from the undiluted extracts as the level showing a signal-to-noise ratio of at least 10 to 1. LOQ for perchlorate was 0.002 mg/kg in cucumber, apple, barley, orange and tomato matrix. The validation in apple and barely matrix (n=5) showed satisfying recoveries between 70–120 % and relative standard deviations (RSDs) well below 10 %. For apple, the average recovery was 108 % at 0.01 mg/kg and 98 % at 0.1 mg/kg. RSDs were 2.6 % and 2.2 %, respectively. For barley, the average recovery was 106 % for the 0.01-level and 93 % for the 0.1-level and RSDs were 2.3 % and 5.0 %, respectively.

The method was applied to 350 conventionally and 78 organically grown samples from the local market. Perchlorate was detected in 66 (19 %) of the conventional samples and in two (2.6 %) of the organic samples. In six of all analyzed samples the PMTDI of 0.01 mg/kg bodyweight was exceeded. The average perchlorate amount of all 66 samples containing perchlorate was 0.043 mg/kg, the highest perchlorate amount detected was 0.4 mg/kg in a tomato sample from Spain.¹

Figure 4 a) and b) show an overview of perchlorate findings in all analyzed fruit and vegetable samples. It was remarkable, that perchlorate could most frequently be detected in melon (50 % of all analyzed melon samples), sweet pepper (53 %), tomato (40 %) and citrus (23 %) commodities. However, the number of analyzed samples is yet not representative enough and does not allow general statements, especially since not all commodity groups have been implemented to a significant number.

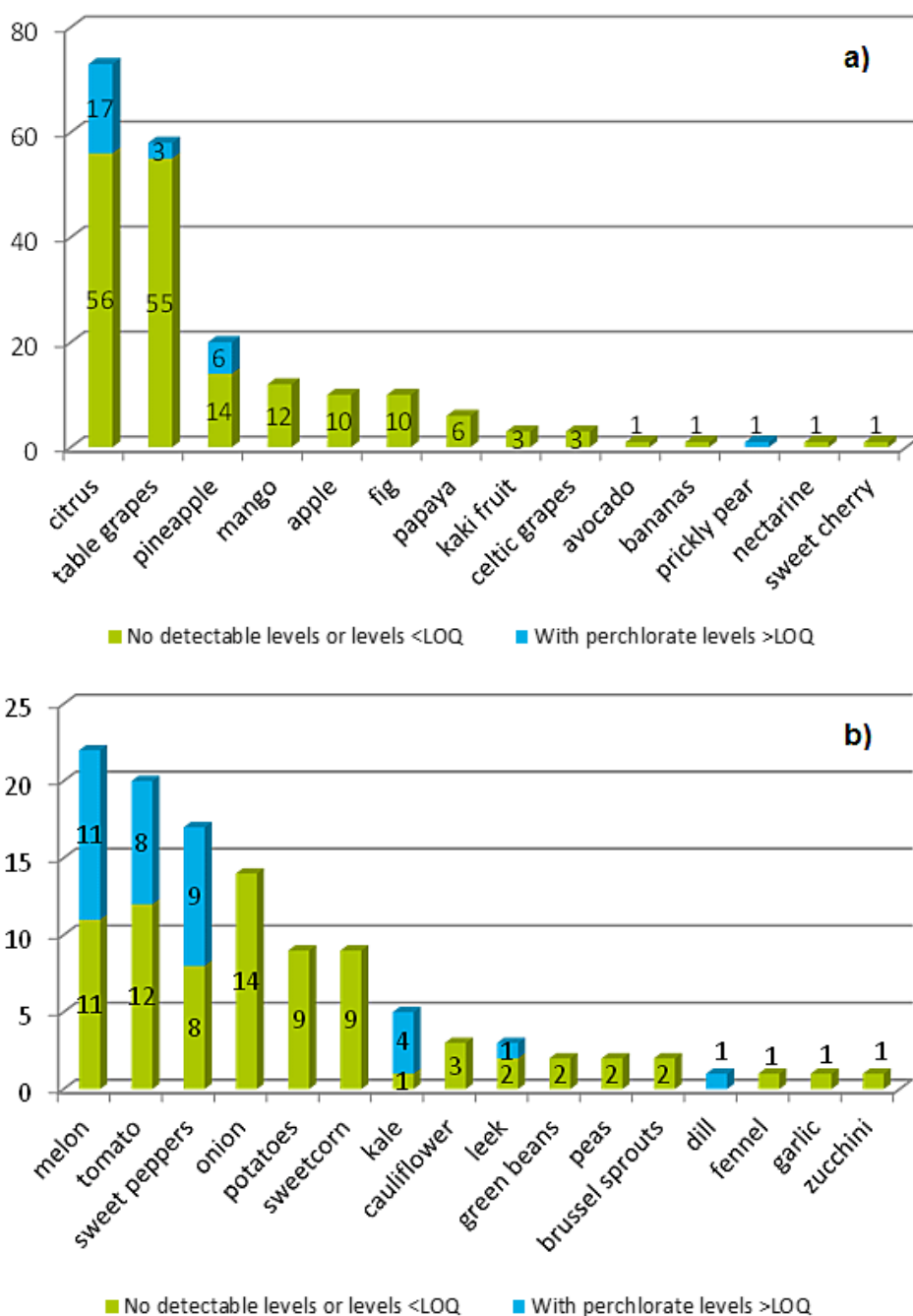


Fig-4: Number of samples with and without perchlorate findings in a) fruit commodities b) vegetable commodities

Figure 5 shows how many samples from each country of origin contained perchlorate, irrespective if of conventional or organic cultivation. 37 % of all samples from Spain and 28 % of all samples originating from Turkey contained perchlorate. In contrast, only 7 % of all samples from Germany contained perchlorate. Once again, the number of analyzed samples does not allow general statements and more samples have to be analyzed for perchlorate.

Table 1 shows in which types of samples perchlorate was found and from which countries these samples originated from. The

highest perchlorate findings detected were 0.4 mg/kg in a tomato sample from Spain, 0.18 mg/kg in a grapefruit from Spain and 0.31 mg/kg in a tomato sample with no specified country of origin.

Conclusions

Our study demonstrates that the QuPPE method is well applicable for the monitoring of the environmental contaminant perchlorate. Determinative analysis is achieved by LC-MS/MS along with a number of highly polar pesticides. Isotope labeled perchlorate is used as internal standard. The

¹The data only includes samples that were analyzed after the method validation was completed.

validation of the method showed satisfying recoveries and RSDs in apple and barley matrix. We recommend the use of the ^{18}O -labeled internal standard, if available. Otherwise the use of matrix matched calibration and the dilution of QuPPE extracts should be mandatory.

Validation experiments for perchlorate additional commodities such as citrus fruit and dry fruit are in process. Moreover, the routine monitoring of perchlorate in fruit and vegetable samples from the local market will be continued.

Acknowledgement:

A special thank you to Mr. Jean-Luc Deborde from SCL de Strasbourg, who provided the introductory information about the perchlorate issue in France.

Table 1: Perchlorate findings in different commodities and their country of origin

Country of origin	Commodity	Perchlorate levels in mg/kg
Spain	tomato (5x)	0.003 – 0.4
	clementine (3x)	0.018 – 0.028
	grapefruit (3x)	0.008 – 0.18
	lemon (2x)	0.011 – 0.017
	mandarin (2x)	0.019 – 0.027
	orange (2x)	0.011 – 0.077
	table wine grape white	0.035
	Turkey	sweet pepper (8x)
melon		0.023
grapefruit		0.014
Brazil	melon (8x)	0.009 – 0.037
Germany	kale (5x)	0.013 – 0.032
	porree	0.059
	table wine grape white	0.012
Peru	clementine	0.064
	mandarin	0.063
	table wine grape red	0.015
Morocco	sweet pepper	0.090
	tomato (2x)	0.010 – 0.043
Costa Rica	pineapple (3x)	0.003 – 0.073
Ghana	pineapple (3x)	0.020 – 0.054
Belgium	kale (2x)	0.019 – 0.022
Chile	mandarin	0.20
Netherlands	frozen pre-fried French fries	0.027
Italy	prickly pear	0.017
China	pomelo	0.010
Honduras	melon	0.010
unknown	dill	0.016
	frozen pre-fried French fries	0.012
	melon	0.015
	tomato	0.31

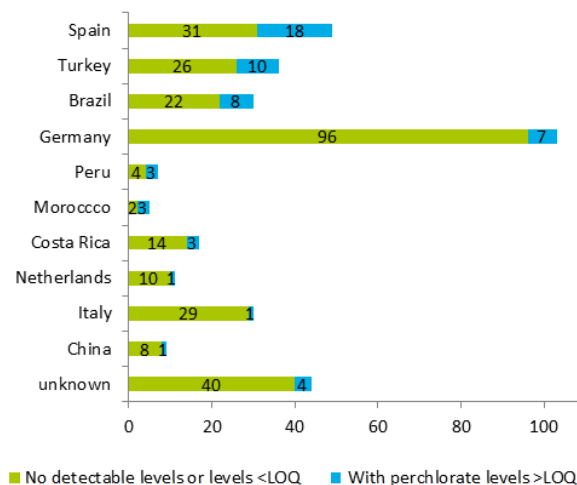


Fig 5: Country of origin - Number of samples with and without perchlorate residues

References

1. B. E. Jander G., *Lehrbuch der analytischen und präparativen anorganischen Chemie (1989)* 157 f., S. Hirzel Verlag, Stuttgart
2. The German Federal Environment Agency (Umweltbundesamt), *Bericht des Umweltbundesamtes vom 18.09.2012 über das Vorkommen und die Verwendung von Perchloraten und deren Vorläufern sowie wesentliche Eintragspfade von Perchloraten in Lebensmittel (2012)* Dessau-Roßlau
3. E. T. Urbansky, *Perchlorate as an environmental contaminant, Environ Sci Pollut Res Int (2002)* 9(3): 187–92
4. COMMISSION DECISION of 10 November 2008 concerning the non-inclusion of chlorate in Annex I to Council Directive 91/414/EEC and the withdrawal of authorisations for plant protection products containing that substance (2008) *Official Journal of the European Union, Brussels*
5. EU Pesticide Database, European Commission Directorate-General Health & Consumers, http://ec.europa.eu/sanco_pesticides Cited 8 April 2013
6. REGULATION (EC) NO 396/2005 OF THE EUROPEAN PARLIAMENT AND OF THE COUNCIL on maximum residue levels of pesticides in or on food and feed of plant and animal origin and amending Council Directive 91/414/EEC of 23 February 2005 (2005)

- Official Journal of the European Union, Strasbourg
7. A. Srinivasan and T. Viraraghavan, Perchlorate: health effects and technologies for its removal from water resources, *Int J Environ Res Public Health* (2009) 6(4): 1418–42
 8. M. A. Greer, G. Goodman, R. C. Pleus and S. E. Greer, Health effects assessment for environmental perchlorate contamination: the dose response for inhibition of thyroidal radioiodine uptake in humans, *Environ Health Perspect* (2002) 110(9): 927–37
 9. Joint FAO/WHO Expert Committee on Food Additives (JECFA), Evaluation of certain contaminations in food (2011) Seventy-Second Report of the Joint FAO/WHO Expert Committee on Food Additives, 64–73, 101
 10. National Academy of Sciences (NAS), Health implications of perchlorate ingestion. Committee to assess the health implications of perchlorate ingestion (2005) The National Academic Press, Washington, DC.
 11. U.S. Environmental Protection Agency, Ion chromatographic determination of perchlorate: analysis of fertilizers and related materials (2001) i-vii, 1–20
 12. R. A. Niemann, A. J. Krynetsky and D. A. Nortrup, Ion Chromatographic Determination of Perchlorate in Foods by On-Line Enrichment and Suppressed Conductivity Detection, *J. Agric. Food Chem.* (2006) 54(4): 1137–1143
 13. D. K. MacMillan, S. R. Dalton, A. J. Bednar, S. A. Waisner and P. N. Arora, Influence of soil type and extraction conditions on perchlorate analysis by ion chromatography, *Chemosphere* (2007) 67(2): 344–50
 14. S. Uchikawa, Spectrophotometric determination of perchlorate ions by solvent extraction with Crystal Violet, *Bull. Chem. Soc. Jpn.* (1967) 40(4): 798–801
 15. A. Soleymanpour, A. Hanifi and A. H. Kyanfar, Polymeric membrane and solid contact electrodes based on Schiff base complexes of Co(III) for potentiometric determination of perchlorate ions, *Bull. Korean Chem. Soc.* (2008) 29(9): 1774–1780
 16. F. Gholamian, M. A. Sheikh-Mohseni and M. Salavati-Niasari, Highly selective determination of perchlorate by a novel potentiometric sensor based on a synthesized complex of copper, *Mater. Sci. Eng., C* (2011) 31(8): 1688–1691
 17. J. Gertsch and C. Henry, Determination of Perchlorate in Drinking and Surface Water by Microchip Capillary Electrophoresis, (2008)
 18. J. C. Gertsch, S. D. Noblitt, D. M. Cropek and C. S. Henry, Rapid analysis of perchlorate in drinking water at parts per billion levels using microchip electrophoresis, *Anal Chem* (2010) 82(9): 3426–9
 19. S. Chattaraj, K. De and A. K. Das, Indirect determination of perchlorate by atomic absorption spectrometry, *Mikrochim. Acta* (1992) 106(3–6): 183–90
 20. W. J. Collinson and D. F. Boltz, Indirect spectrophotometric and atomic absorption spectrometric methods for determination of perchlorate, *Anal. Chem.* (1968) 40(12): 1896–8
 21. T. W. Collette and T. L. Williams, Determination of perchlorate in some fertilizers and plant tissue by Raman spectroscopy, (2001)
 22. M. Kamoshita, K. Kosaka, M. Asami and Y. Matsuoka, Analytical method for perchlorate in water by liquid chromatography-mass spectrometry using an ion exchange column, *Anal Sci* (2009) 25(3): 453–6
 23. E. Hedrick and D. Munch, Measurement of perchlorate in water by use of an 18O-enriched isotopic standard and ion chromatography with mass spectrometric detection, *J Chromatogr A* (2004) 1039(1-2): 83–8
 24. J. Romano, The determination of perchlorate in water using LC-MS-MS, *Am. Lab.* (Shelton, CT, U. S.) (2006) 38(17): 12–14
 25. P. Winkler, M. Minter and J. Willey, Analysis of perchlorate in water and soil by electrospray LC/MS/MS, *Anal Chem* (2004) 76(2): 469–73
 26. E. Burneson, Drinking water: regulatory determination on perchlorate, *Fed. Regist.* (2011) 76(29): 7762–7767
 27. S. M. Backus, P. Klawuun, S. Brown, I. D'sa, S. Sharp, C. Surette and D. J. Williams, Determination of perchlorate in selected surface waters in the Great Lakes Basin by HPLC/MS/MS, *Chemosphere* (2005) 61(6): 834–843
 28. H.-C. Chen, W.-T. Chen and W.-H. Ding, Determination of perchlorate in river by ion-pair hollow-fiber liquid-phase microextraction coupled with electrospray ionization tandem mass spectrometry, *Talanta* (2009) 79(2): 442–445
 29. R. A. Niemann and A. J. Krynetsky, Analytical methods for the determination of perchlorate anion in foods, (2006)
 30. A. J. Krynetsky, R. A. Niemann and D. A. Nortrup, Determination of Perchlorate Anion in Foods by Ion Chromatography-Tandem Mass Spectrometry, *Anal. Chem.* (2004) 76(18): 5518–5522
 31. A. J. Krynetsky, R. A. Niemann, A. D. Williams and M. L. Hopper, Rapid determination of perchlorate anion in foods by ion chromatography-tandem mass spectrometry, *Prepr. Ext. Abstr. ACS Natl. Meet., Am. Chem. Soc., Div. Environ. Chem.* (2005) 45(2): 914–919
 32. M. Anastassiades, D. I. Kolberg, D. Mack, C. Wildgrube, I. Sigalova, D. Roux, D. Fügler EU Reference Laboratories for Residues on Pesticides - EURL for Single Residue Methods, Quick Method for the Analysis of Residues of numerous Highly Polar Pesticides in Foods of Plant Origin involving Simultaneous Extraction with Methanol and LC-MS/MS Determination (QuPPE-Method) - Version 7 (December 2012) http://www.crl-pesticides.eu/library/docs/srm/meth_QuPPE.pdf Cited 3 April 2013
 33. Datapool of the EU Reference Laboratories for Residues of Pesticides, www.eurl-pesticides-datapool.eu Cited 10 April 2013
 34. H. Stahnke, S. Kittlaus, G. Kempe and L. Alder, Reduction of matrix effects in liquid chromatography-electrospray ionization-mass spectrometry by dilution of the sample extracts: how much dilution is needed?, *Anal Chem* (2012) 84(3): 1474–82